

# Single-Use Urology Videoscope

MDA Reg No. : GB10740222-103413

**Working Channel:**  
3.6Fr (Φ1.2mm)



## More Precise

The maximum 60° rotation function allows the operator to change the shafts angle by rotating. This enables more user precision and comfort.


## Field of View: 120°

Wider field of view of 120° gives the operator a better visual experience



## Deflection 285° UP/DOWN

Deflection capabilities of 285° up and down allow for optimal visualization throughout the whole urological tract.

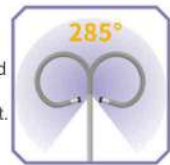

## The uniform outer diameter 8.4Fr(Φ2.8mm)

Slim, pliable and tough insertion core to help operator navigate smoothly in urethra and kidney, ensuring sufficient irrigation flow.

## The rotation function

At the top of the device allows the operator to manage the operation in a standing or sitting position

## Excellent passive bending

Effectively reach the lower calyces



**No sterilization required**



**No cross-infection risk**




**No maintenance required**

Three different type of irrigation valve satisfy different using habit, ensuring a safe fit for laser fibers and reducing leakages.

80460



T-001

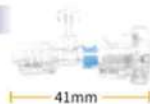

Y-001



## Be ready at any time

- Plug and play. Off the Shelf Ready. Consistency every time.
- You decide which hand drives. Avoid cross-infection.
- Save on re-sterilisation and maintenance costs.
- Eradicate unnecessary repair bills.
- Save time post procedure, no need for scope integrity testing.
- Angulation lock function
- Unlimited operational Time



**Authorised Representative**  
IBS MEDICAL EQUIPMENTS SDN BHD  
No. 7-G, 3rd Mile Square Business Centre,  
No. 151, Old Klang Road, 58100 Wilayah  
Persekutuan Kuala Lumpur.



Tel : +603 7971 4221



Whatsapp : +6012 346 3400




Website : [www.ibsmedic.com](http://www.ibsmedic.com)




Email : [sales@ibsmedic.com](mailto:sales@ibsmedic.com)

# Endoscopic Image Processor

MDA Reg No. : GB6207922-110588

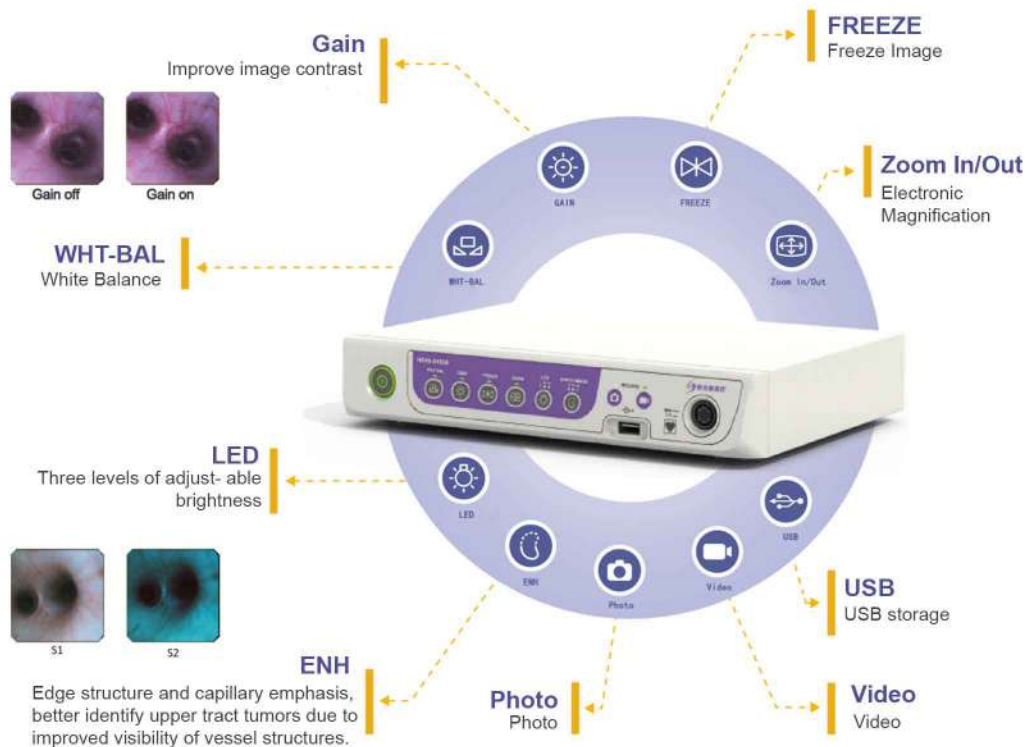

Scivita brings a brand-new and higher-resolution digital image without distortion, even on a 32-inch monitor.

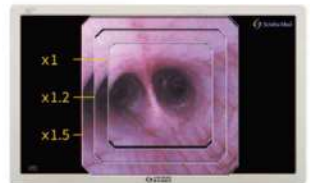

The processor compatible with various monitors. 32-inch monitor is highly recommended.

## Convenient for doctors

- Hold horizontally or vertically, easy to operate
- Ergonomic design

## Beneficial for patients

- Plug and play, increase turnover rate
- No maintenance required

## Beneficial for patients

- Plug and play, increase turnover rate
- No maintenance required

## SPECIFICATION



	Single-Use Ureteroscope		Single-Use Cystoscope
Model	SUV-2A-B	SUV-2C-B	SUV-1D-B
Shaft Diameter OD	8.4Fr ( $\Phi$ 2.8mm )	7.5 Fr ( $\Phi$ 2.5mm)	16.2Fr ( $\Phi$ 5.4mm)
Outer Diameter of distal End	7.8Fr ( $\Phi$ 2.4mm )	7.2 Fr ( $\Phi$ 2.4mm)	15.6Fr ( $\Phi$ 5.2mm)
Weight	200g	200g	200g
Working Length	680mm	680mm	380mm
Total Length	990mm	990mm	700mm
Working Channel OD	3.6Fr ( $\Phi$ 1.2mm )	3.6 Fr ( $\Phi$ 1.2mm)	6.6Fr ( $\Phi$ 2.2mm)
Field of View	120°	120°	120°
Depth of Field	0.5~120mm	0.5~120mm	0.5_120mm
Angulation Range	285° Up/Down	285° Up/Down	220° Up/ Down

## SPECIFICATION

### PROCESSOR



Model	HDVS-S100A
Size	300mm X 225mm X 57mm
Weight	2.5 Kg
Power Supply Requirement	100-240V
Output	DV*1, SDI*2
Interface	Foot Switch and USB

### ACCESSORIES(OPTIONAL):

Irrigation Valve	80460	39mm
	T-001	41mm
	Y-001	84mm
Foot Pedal	KF2 1S-MED GP25	

